


Liquid Dividers

The manually operated column heads are used to adjust the reflux ratio between reflux and drain in distillation columns. Depending on the version, they are designed for operation under vacuum and/or normal pressure.

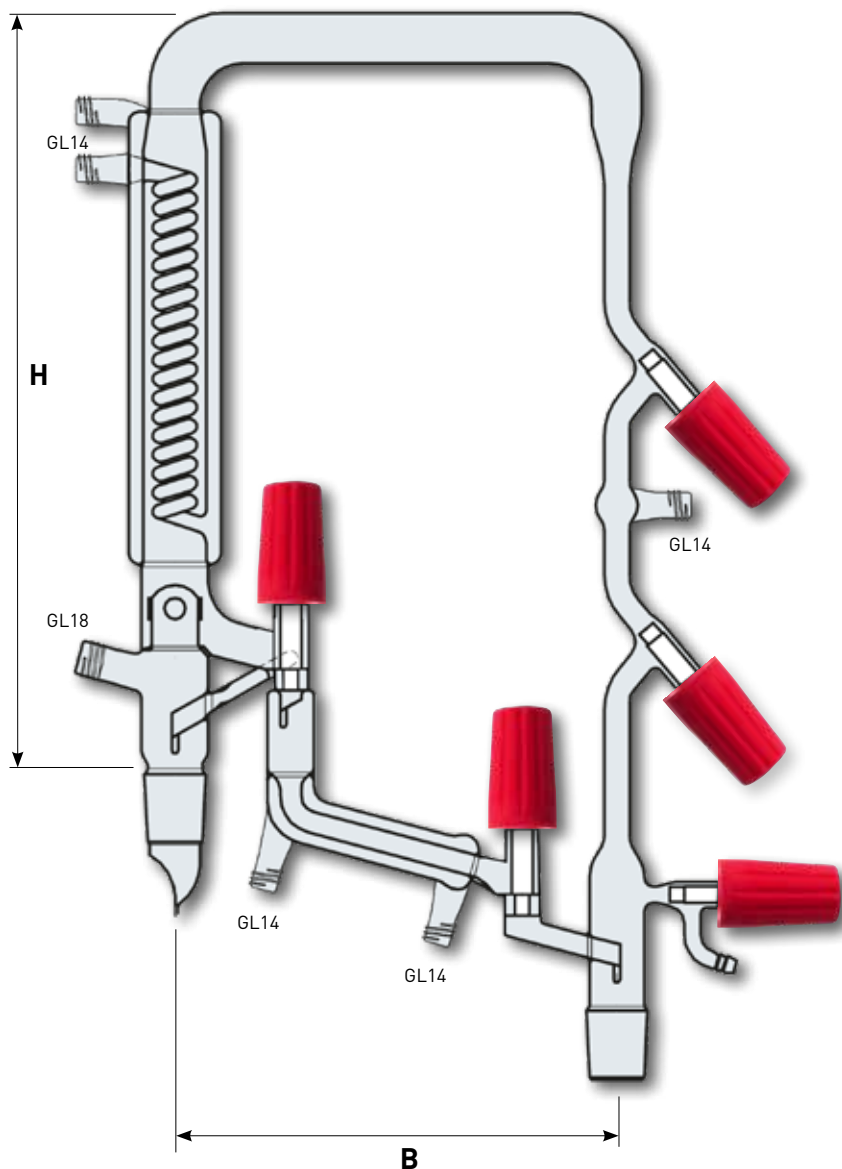

As there are countless variants and variations, we present here our best-selling liquid divider. Of course, all parameters of this model can also be configured according to your wishes. We will gladly offer you a liquid divider according to your needs.

PRODURAN® Valves



- Precise dosing
- Unscrew protection
- Vacuum-tight, leak rate $< 1 \times 10^{-6}$ mbar x l/s

Liquid Divider with Intensive and Intercooler, Dead Volume free

In- an outlet NS 29/32, intensive cooler 200 mm

Article No.	H mm	B mm	Valves	Vakuum operation
20.738.30	370	170	Produran 6	yes



All GL-14 connections with cap and PP olive (optionally PTFE)